


1 Installing EzDent-i




A **STEP 1:** Select installation file "Setup.exe"

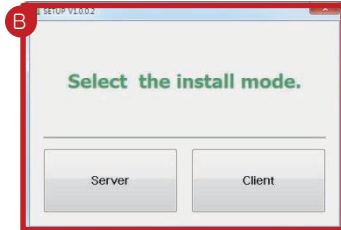

B **STEP 2:** The install mode selection window will appear

Server: Installs the file server(EzServer), creates the patient database (DB), and the imaging software(EzDent-i)

Client: Installs only the imaging software (EzDent-i)

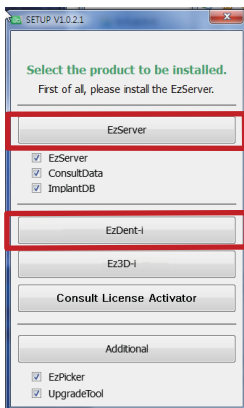

C **STEP 3:** The install mode selection window will appear

EzServer: Installs the file server, and creates patient database (DB)

ConsultData: Installs VCT(Virtual Consultation Tool)

ImplantDB: Installs implant fixture data for implant simulation

D **STEP 4:** Installs the imaging software (EzDent-i)

! **NOTE:** Console Vs Client version of EzDent-i



Console Version


This version of EzDent-i has the Extraoral X-ray device(CBCT, Pano, Ceph) capture function enabled, this is meant for capture PC installations

Client Version


This version of EzDent-i **DOES NOT** have the Extraoral X-ray device (CBCT, Pano, Ceph) capture function enabled, this is meant for Intraoral X-ray device installations

License Activation With A Serial Key (Softkey Method)


2

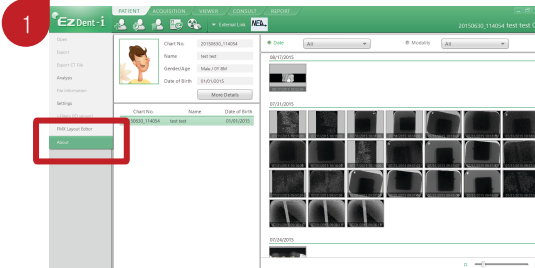


Start EzDent-i, click on EzDent-i logo on the top-left corner, click [About] from the dropdown menu

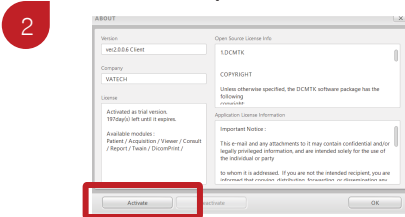

Click [Activate] and begin license activation

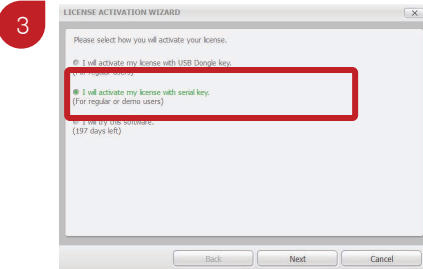

Select "I will activate my license with serial Key" and click [Next]

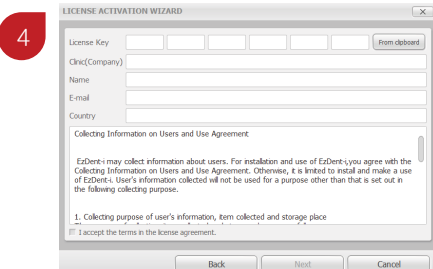


Insert the serial license key and click [Next], activation is complete when an "Activation Successful" prompt appears


License Activation With A USB Dongle Key (VTNL Method)




Insert the USB Dongle key into the PC, locate and execute the "VTNL.exe" file contained within the Dongle Key

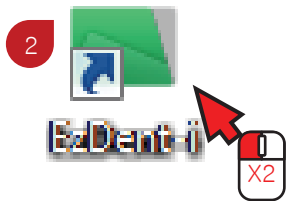

Click and start EzDent-i

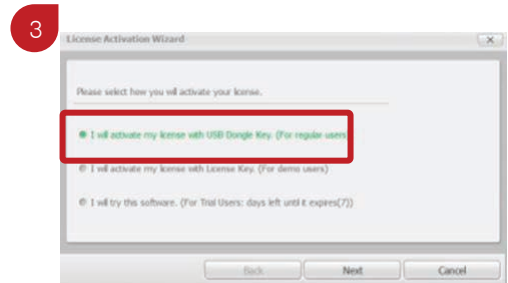

Select "I will activate my license with USB Dongle Key" and click [Next]

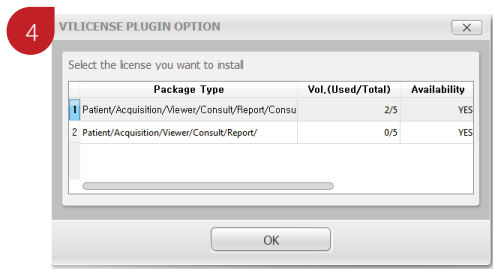

Select the license key from the list and click [OK], activation is complete when an "Activation Successful" prompt appears